



Microscope Cleaning



- Blow off dust with canned air

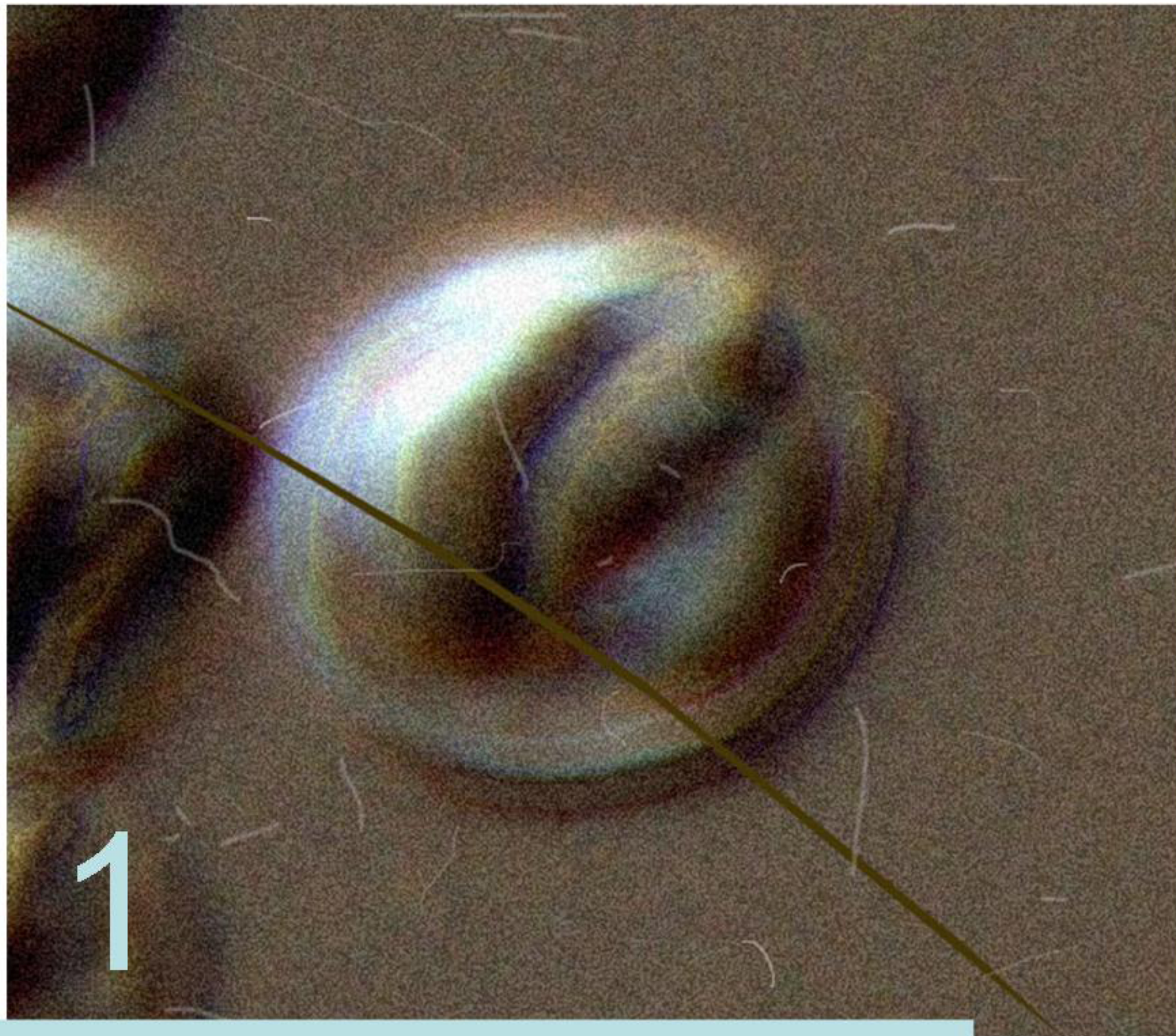

- Clean the oculars
 - Use new, dust-free cotton-tipped stick dampened with microscope cleaning solution
 - Start in the center of the lens and spiral to the outside while rotating the cotton tip
 - Repeat with new sticks until the view is clear

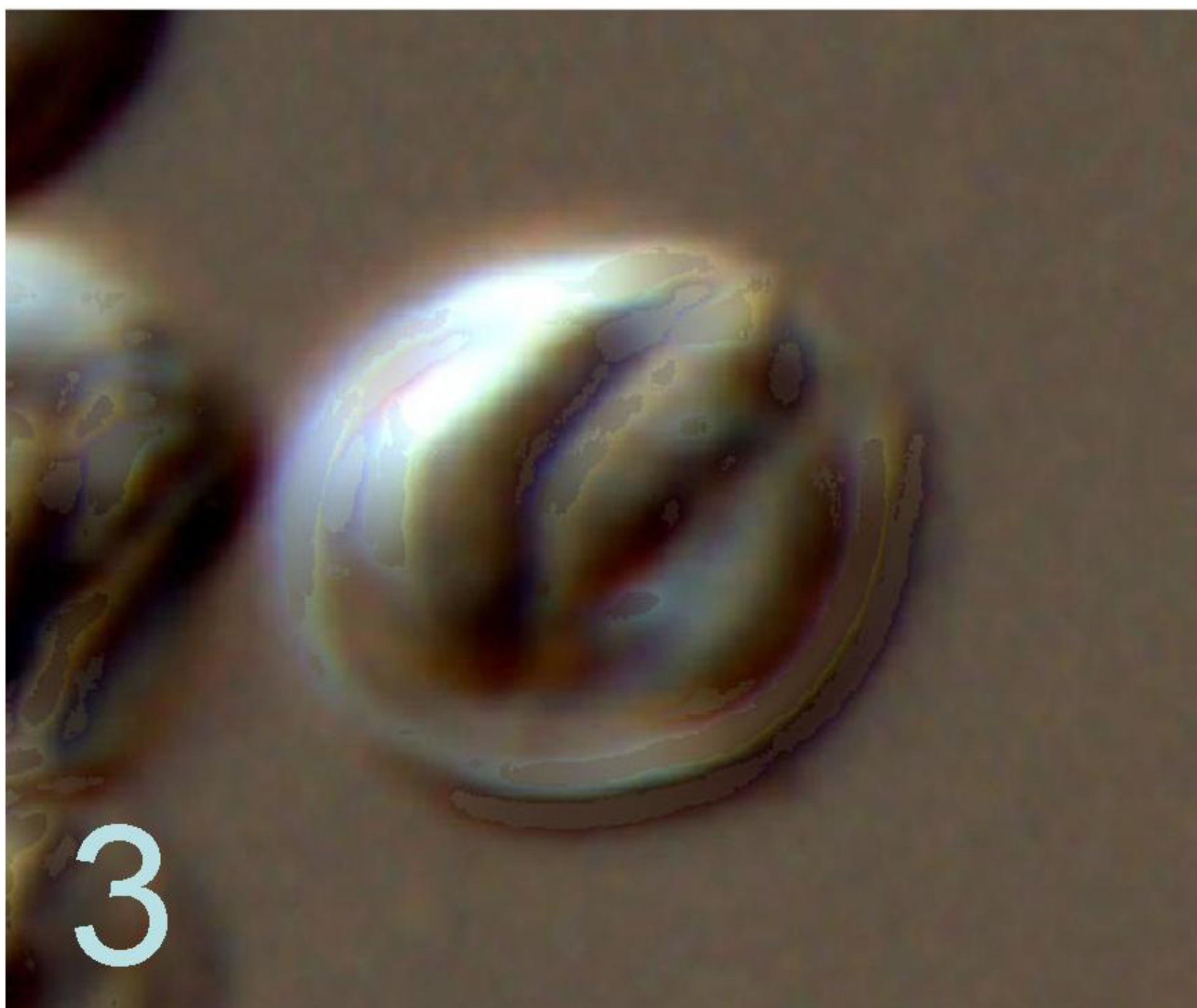

- Clean all the objectives
 - Wipe 100X oil objective with lens paper
 - Hold a new, dust-free cotton-tipped stick dampened with microscope cleaning solution at a 45° angle on the objective lens and twirl
 - Repeat with new sticks until the view is clear



- Repeat steps 1-3 after each session
- Reminder:**
- Thoroughly clean all surfaces and additional lenses weekly
 - Complete annual service per manufacturer guidelines